


*LeTube*TM

Single-use Biliary Draining Tubing

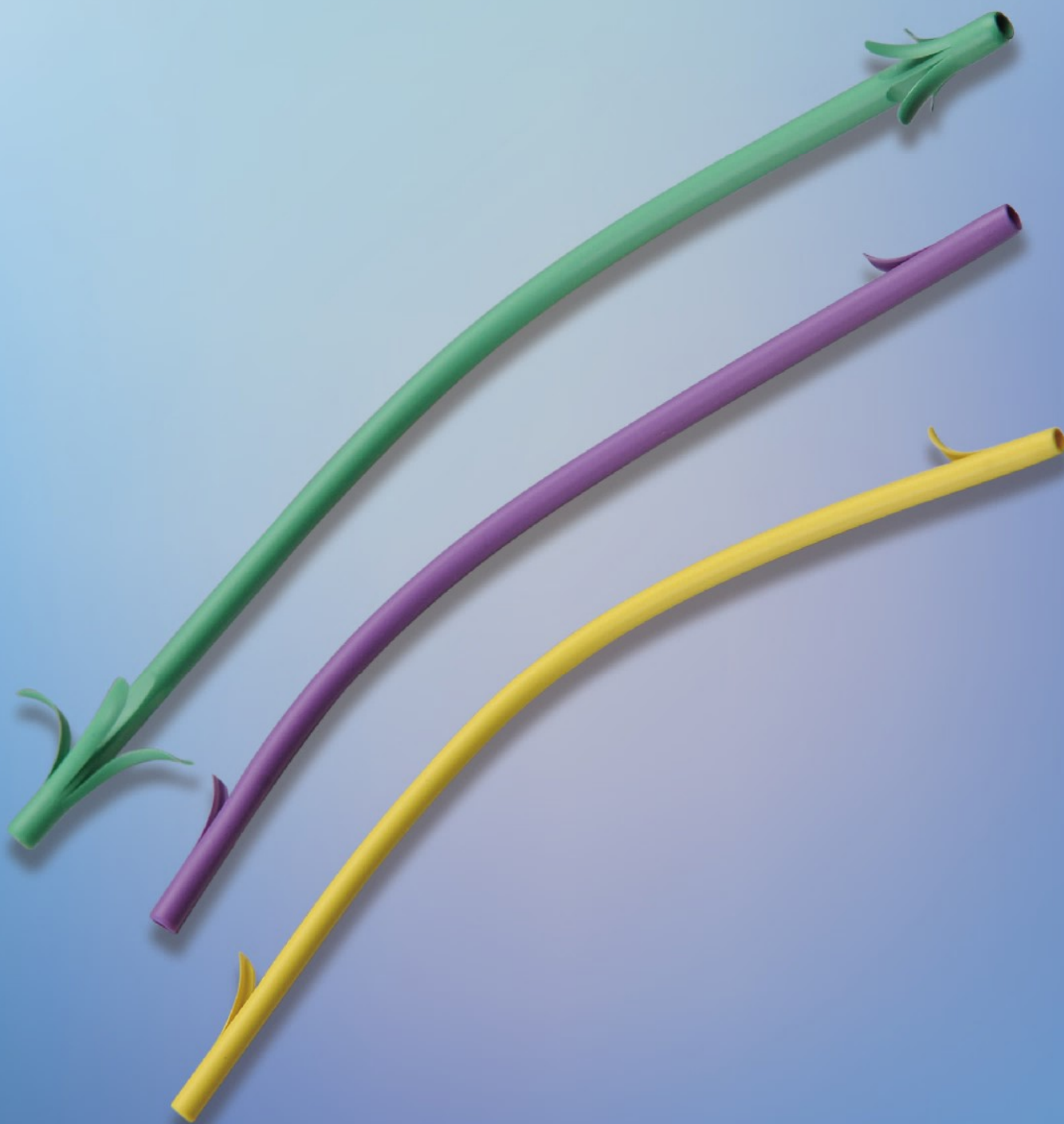

Leo+Med
乐奥医疗 为生命健康服务

Endoscopic Accessories ERCP Series 


LeTube™ Single-use Biliary Draining Tubing

Used for endoscopic retrograde biliary drainage(ERBD).



● Good Radiopacity

- Whole drainage tube has clear contrast image under X-ray during the procedure.
- Multiple markers of delivery system provide enough radiopacity to facilitate position confirmation.

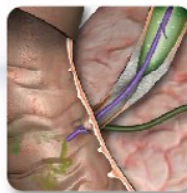

● Precision placement

- Perfect balance of drainage tube flexibility and stiffness improves good deliverability.
- The flaps are designed to increase reliability and reduce migration.
- Preloaded type is designed to provide endoscopic adjustable and precise placement.

● Excellent surface status

- Fluoroplastic material provides lower inner surface friction coefficient compared to conventional plastic materials.
- Excellent surface provides less possibility of bile duct obstruction.

Specification:

▲ Single Wing

Tubing Single REF	Pre-loaded REF	Tubing Length (mm)	Tubing Diam. (Fr)	Tubing Configuration	Compatible Guidewire	Minimum Working Channel (mm)
WL-7-700	WL-7-700-P7	70	7.0	Straight	0.035"	2.8
WD-7-900	WD-7-900-P7	90		Proximal bend	0.035"	2.8
WC-7-1200	WC-7-1200-P7	120		Center bend	0.035"	2.8
WL-8.5-700	WL-8.5-700-P8.5	70	8.5	Straight	0.035"	3.2
WD-8.5-900	WD-8.5-900-P8.5	90		Proximal bend	0.035"	3.2
WC-8.5-1200	WC-8.5-1200-P8.5	120		Center bend	0.035"	3.2
WL-10-700	WL-10-700-P10	70	10	Straight	0.035"	3.8
WD-10-900	WD-10-900-P10	90		Proximal bend	0.035"	3.8
WC-10-1200	WC-10-1200-P10	120		Center bend	0.035"	3.8

▲ Not all part numbers shown on this product information sheet. Please consult with LEOMED representative for details.

▲ Christmas tree

Tubing Single REF	Pre-loaded REF	Tubing Length (mm)	Tubing Diam. (Fr)	Tubing Configuration	Compatible Guidewire	Minimum Working Channel (mm)
CL-8.5-700	CL-8.5-700-P8.5	70	8.5	Straight	0.035"	3.2
CD-8.5-900	CD-8.5-900-P8.5	90		Proximal bend	0.035"	3.2
CC-8.5-1200	CC-8.5-1200-P8.5	120		Center bend	0.035"	3.2
CL-10-700	CL-10-700-P10	70	10	Straight	0.035"	3.8
CD-10-900	CD-10-900-P10	90		Proximal bend	0.035"	3.8
CC-10-1200	CC-10-1200-P10	120		Center bend	0.035"	3.8

▲ Not all part numbers shown on this product information sheet. Please consult with LEOMED representative for details.